

# 3D KNEADING MASSAGER + HEAT



*Acumag*

Operation Instructions  
Read instructions before using

## TABLE OF CONTENTS

Warning and Caution.....	1-4
Control unit .....	4
Using the massager.....	5
Product Care and Maintenance.....	6

## WARNING

### FOR ADULT USE ONLY

Important: Any individual who may be pregnant, has a pacemaker, suffers from phlebitis and/or thrombosis, is at an increased risk of developing blood clots, or who has pins/screws/artificial joints or other medical devices implanted in his/her body should consult with a physician before using a massaging device designed for home use. If HEAT is too strong for your comfort please switch the heat off.

- Do not use on an infant or invalid or on a sleeping or unconscious person.
- Do not use on sensitive skin or on a person with poor blood circulation.

## CAUTION

- TO REDUCE THE RISK OF ELECTRIC SHOCK, DO NOT REMOVE COVER. THERE ARE NO SERVICEABLE PARTS INSIDE.
- TO REDUCE THE RISK OF FIRE OR ELECTRIC SHOCK, DO NOT EXPOSE THIS UNIT TO RAIN OR MOISTURE

## **IMPORTANT SAFETY INSTRUCTIONS**

All of the safety and operating instructions should be read, adhered to and followed before the unit is operated.

When using an electrical appliance, basic precautions should always be followed, including the following:

**DANGER** To reduce the risk of electric shock:

1. Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
2. Do not reach for an appliance that has fallen into water. Unplug immediately.
3. Do not use while bathing or in a shower.
4. Do not place or store appliance where it can fall or into a tub or sink.
5. Do not place in or drop into water or other liquid.
6. Never use pins or other metallic fasteners with this appliance.
7. Carefully examine the covering before each use. Discard the appliance if the covering shows any signs of deterioration, such as checking, or cracking.
8. Keep dry. Do not operate in a wet or moist environment.

**WARNING** - To reduce the risk of burns, fire, electric shock, or injury to persons:


1. An appliance should never be left unattended when plugged in. Unplug from the outlet when not in use, and before putting on or taking off parts.
2. Close supervision is necessary when this appliance is used by, on or near children, invalids, or disabled persons.
3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this appliance if it has a damaged cord or plug. If it is not working properly, if it has been dropped or damaged, or dropped into water, DO NOT attempt to fix it yourself. Return the appliance to a service center for examination and repair.
5. Do not carry this appliance by its supply cord or use the cord as a handle.
6. Do not crush or fold this appliance when storing.
7. Keep the cord away from heated surfaces.
8. Never drop or insert any object into any opening.
9. Do not use outdoors. This appliance is intended for household and indoor use only.
10. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
11. Do not operate in the presence of explosive and/or flammable fumes.

12. To disconnect, set all controls to the off position, then remove plug from the outlet.
13. Do not overload the electrical outlet. Use only the power source as indicated.
14. To avoid the risk of electric shock, do not disassemble or attempt to repair the appliance. Incorrect repair can cause risk of electric shock or injury to persons when the unit is used.
15. Never remove the plug from the outlet by pulling the power cord.
16. Do not use this product as a head tapper.

## SAVE THESE INSTRUCTIONS

### LOCATION OF PARTS AND CONTROLS

1. Main Unit
2. Shiatsu Nodes
3. Control Panel
4. power on/off
5. Massage direction
6. Speed Control
7. Heat on/off
8. Adapter
9. Car Power Plug

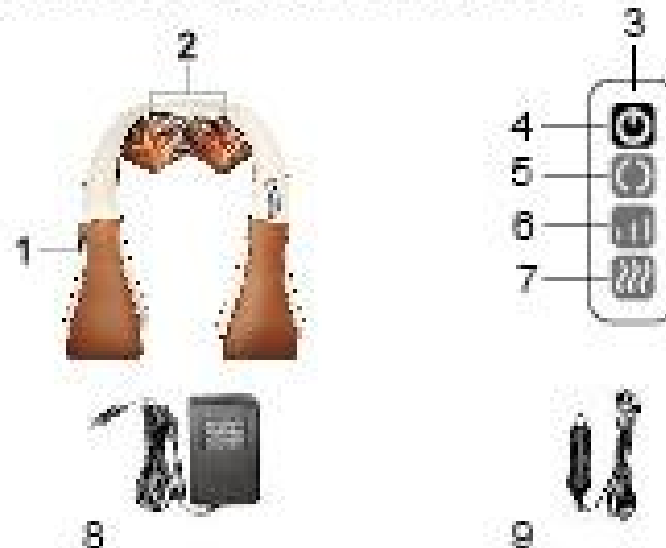

## **USING THE MASSAGER**

1. For **INDOOR USE**, connect the power adaptor to the massager. Plug the adaptor into an outlet. (Sit comfortably in the chair then use massager)
2. For **IN CAR USE**, connect the car power plug to the massager. Plug the car power plug into the cigarette lighter in the car. (Do not use while you are driving)
3. Grip the **POSITIONING HANDLES** and drape the **MAIN UNIT** around your neck with the **MESSAGE NODES** in the center of your back.
4. Press the **POWER** button on/off to turn on/off massage. Press the **HEAT** button on/off to turn on/off heat
5. Press the **HEAT** button on/off to turn on/off heat.
6. Press **DIRECTION** button to change the direction of massage nodes, clockwise or counterclockwise.
7. Press speed control to change the speed

### **NOTE:**

- The massage function will automatically shut off after 10 minutes.
- The Massage direction will automatically change every one minute.



## **PRODUCT CARE AND MAINTENANCE**

1. Place the massager in a safe; cool and dry place when not in use.
2. After unplugging, wipe with a soft, damp cloth using non-abrasive cleaners.
3. Never submerge the unit in any liquid.
4. Keep away from all solvents and harsh detergents.
5. **DO NOT** attempt to repair this massager. There are no user-serviceable parts.



**ACUMAG LTD**  
**UNIT G1, EUROPA HOUSE**  
**FRAISER ROAD, ERITH**  
**KENT DA8 1QL, UK**  
**[www.acumag.co.uk](http://www.acumag.co.uk)**  
**Tel. 0044 1322 447610**

## **3D KNEADING MASSAGER + HEAT**

### **SPECIFICATIONS**

Rated Voltage.....	DC12V
Frequency.....	50-60Hz
Power.....	24W